

AeroChamber* MV Holding Chamber

AeroChamber* MV Holding Chamber is designed for delivery of metered dose inhaler (MDI) medications in mechanical ventilation circuits or with a manual resuscitation bag.

Refer to the [instructions for use](#) for additional information.



- Clinically validated for use in ventilator circuits¹

Use with Resuscitation Bag

Use in Mechanical Ventilation


22mm Inner Diameter Connector

Patient Wye

Transparent Chamber
Clear chamber allows for visualization of aerosol

Canister Port
Can be used with most commonly prescribed metered dose inhalers (excluding those with integrated dose counters)

Inspiratory Limb

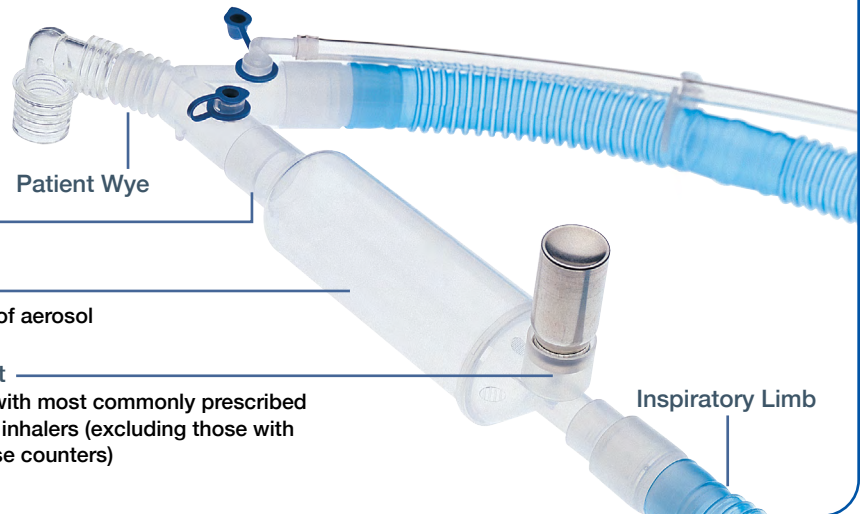

Part Number	Case Quantity
825110050	50

SINGLE-PATIENT – MULTIPLE USE | Replace with the ventilator circuit. Do not clean or disinfect.
PACKAGING: Flow-wrapped with instructions for use.
Not made or manufactured with bisphenol A (BPA), phthalates, latex or lead.

1. [EFFECTS OF INHALED FORMOTEROL COMPARED WITH SALBUTAMOL IN VENTILATED PRETERM INFANTS](#), Rieger-Fackeldey E, et al. Pulmonary Pharmacology & Therapeutics 2004;17(5):293-300.

AeroChamber* MV
Holding Chamber



+1-519-455-4862
info@trudellmed.com



Designed and
Developed in Canada
Manufactured in Canada
with USA parts.

trudellmed.com